


# A choice that matters



# A universal solution

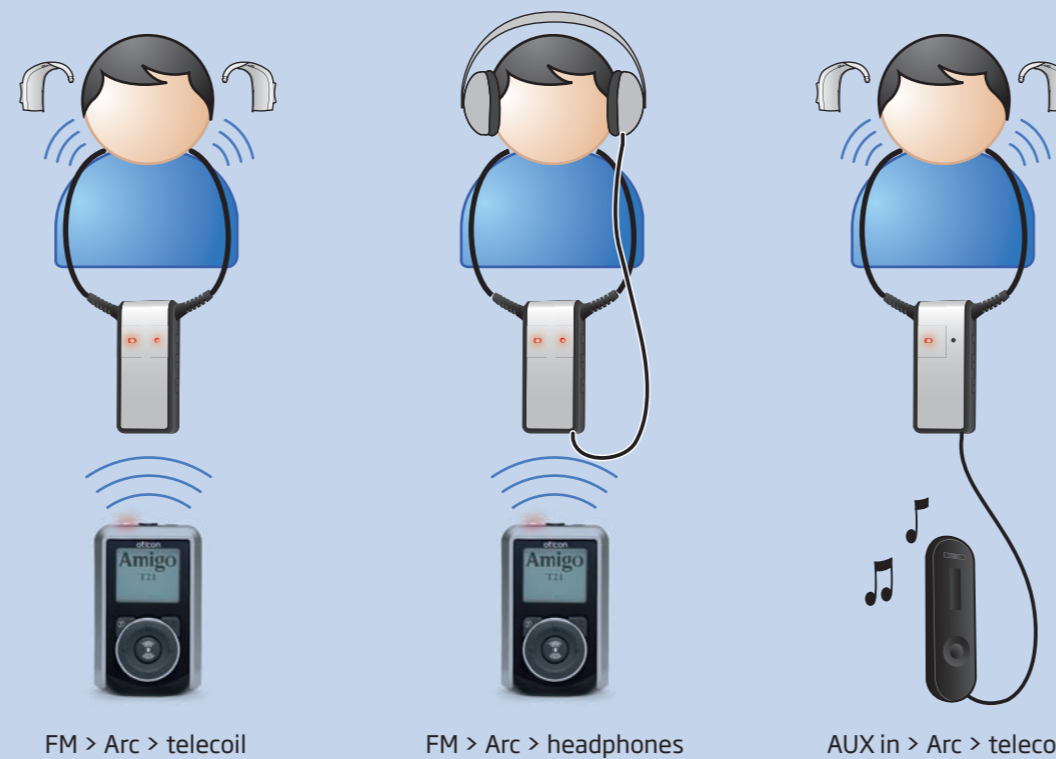

# Features that make a difference

FEATURE	BENEFIT	WHY THIS IS IMPORTANT
7000 Hz bandwidth	Wide bandwidth	Delivers the broadest signal for natural sound quality and enhanced speech intelligibility.
Intelligent sleep mode	Auto sleep mode when not in use	Allows the child to use Amigo Arc all day long.
Battery and FM LEDs	Visual status indicator	Both caregiver and child can be certain that Amigo Arc is working as it should.
USB charging capabilities	Flexible charging options	Easy, instant charging when necessary.
Intelligent volume control	Remembers previous volume settings	Both caregiver and child can trust that the volume settings will remain comfortable.
Works with Amigo T20/T21/WRP that doubles as wireless synchronizer/programmer	Wireless synchronization and programming	Synchronization and programming happens instantly, saving precious time.
Tamper resistant controls	Buttons can be programmed to be active or disabled	Caregivers can be certain the settings are correct.
Low frequency gain emphasis (+6dB)	Superior bass boost	Bass boost in the low frequencies enriches the experience when listening to music.
AUX programmed to IN or OUT	Versatile AUX function	Flexibility to connect either headphones or an external audio source.
Five exchangeable colored shells	Customizable design	Creates a sense of ownership for children, making them feel more motivated to wear Amigo Arc.

### What is Amigo Arc?

Amigo Arc is a universal neckloop receiver that acts as a bridge between an FM transmitter and the hearing aid. Amigo Arc works with any instrument with a telecoil/mic+telecoil, such as BTE, RITE, ITE, cochlear implants or bone anchored instruments.

### Why Amigo Arc?

Although FM has proven to be beneficial in classrooms, connecting these systems to hearing aids is not always easy. Amigo Arc, however, is simple and intuitive to use, offering consistent access to speech, day after day.

### Higher bandwidth - more natural sound

For speech intelligibility, nothing beats Amigo Arc. Its extended, 7000 Hz bandwidth captures far more of the important high frequency speech cues than neckloops with a range of only 5500 Hz.

### More options - higher flexibility

Amigo Arc offers superior sound quality not only for FM but also for music. It has a unique AUX IN option, where the bandwidth increases to 9000 Hz - plus a bass boost that promises far richer sound quality than any other neckloop currently available.

Another exclusive feature is the selection of colorful, exchangeable covers that offer children a chance to put their own personal touch on Amigo Arc.



### Making life easier at school

By acting as a bridge between the FM transmitter and the telecoil-equipped instruments, Amigo Arc offers children a significant advantage in class. The extra boost can also aid concentration in children with normal hearing, e.g. children with auditory processing disorders or attention difficulties.

A very useful feature in Amigo Arc is the visual status indicator (LED), which provides instant confirmation that the FM system is working and the child is hearing clearly. This makes both caregiver and child feel more confident.

### Making life at home more entertaining

Whenever children want to listen at home they can use Amigo Arc in one of two ways: they can connect it directly to a music player or computer, or they can pick up the signal from an FM transmitter that's connected to their TV. This allows them to listen at the same volume as the rest of the family.

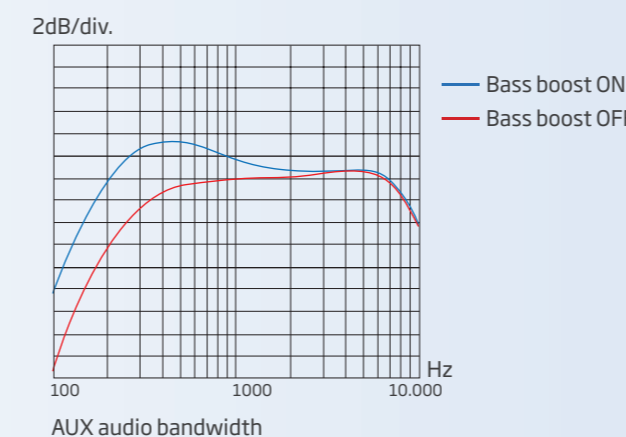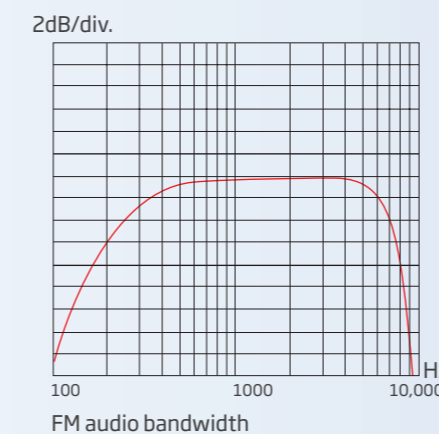

\*applies to Oticon wireless-enabled hearing instruments

# The Amigo family at a glance



**child friendly hearing care**

It takes a truly dedicated approach to help children with hearing problems achieve their full potential. That's why we deliver the solutions and services that professionals and caregivers need to provide children the opportunities they deserve. This is what child-friendly hearing care is all about.



**Amigo Arc**

More of what matters in FM

- Amigo T21/T20:** Advanced, robust FM transmitter and channel synchronizer in one. With Amigo T21 a built-in team teaching function is available.
- Amigo T10:** This elegant, multi-purpose FM transmitter offers three microphone modes: Omnidirectional, Focus and Super Focus.
- Amigo T5:** This sensible, easy-to-use transmitter fulfils all the basic FM requirements.
- Amigo WRP:** Optional accessory for programming and synchronizing Amigo Arc.
- Amigo eZync:** This wall-mounted channel synchronizer automatically synchronizes Amigo Arc and other Amigo receivers as students pass by.



Black Dark red Dark blue White Silver  
Personalize Amigo Arc with colorful shells that are easy to click on and off.