

Amigo FM products are refreshingly easy to use - and the latest addition to this family is no exception. The new Amigo eZync is a dedicated channel synchronization unit that makes channel management easier in educational settings. By automatically synchronizing students' receivers, Amigo eZync saves precious teaching time.

Channel Management made easy

With Amigo eZync installed in classrooms or institutions, Amigo ear-level receivers will be automatically synchronized to the desired channel.

Neat and discreet

Simple to install, Amigo eZync looks neat on a wall or on a shelf. Changing channels is easy using the wireless programming feature of Amigo T20/T21 or WRP.

Instant confirmation


The Amigo eZync is equipped with a unique LED (Light Emitting Diode), which provides instant visual confirmation that Amigo eZync is turned on and synchronization has been successful. The LED also indicates when the Amigo eZync has switched to back-up battery power, and when it has accepted a channel change command.

Amigo eZync promises

- More time for teaching
- Easy installation
- Simple operation
- Small and discreet design
- Maintenance free

Key Features

- LED status indicator
- Built-in battery back-up (up to 48 hours), no need for battery replacement
- Wireless synchronization
- 5 ft., 360 degree synchronization range
- Change of channel is restricted to the first 60 seconds after Amigo eZync has been turned on to prevent accidental changes
- Can be mounted on a wall or stand on a table or shelf

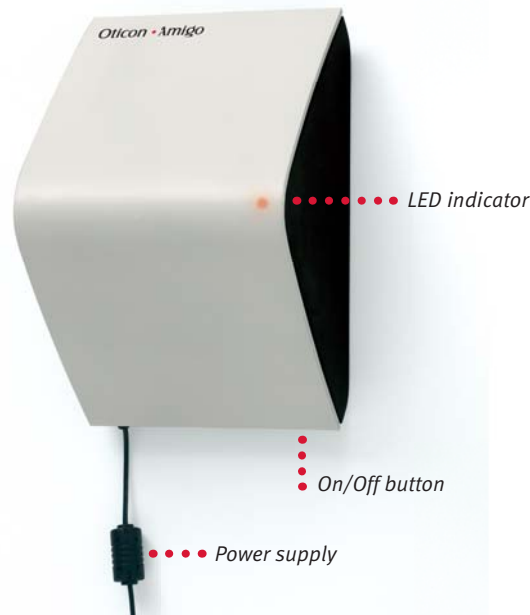

Standard accessories

- Shelf stand
- USB connector AC power supply
- Extension cable for the power supply
- Wall mount
- Fasteners and anchors

Compatibility with other FM systems

Amigo eZync will automatically synchronize Oticon Amigo ear-level receivers.

Technical data

Dimensions (Length x Width x Depth)	19.7 cm x 12.4 cm x 11.1 cm
Weight	1 lb. (459 g)
Power supply	5V, 75 mA
Battery back-up type	NiMH
Battery back-up operation	Up to 48 hours
Charging time	24 hours
Operating frequencies	3.7 MHz
Working range	5 ft. (1.5 m), 360 degree synchronization range
RF output power	Less than -15 dBμA/m @ 10 m
Type of modulation	AID
Antenna type	Coil antenna
Emissions	Amigo eZync complies with all relevant standards